



# OOSAFE® IVF LABORATORY DISINFECTANTS

used by leading IVF clinics in more than 100 countries



GOosafe in IVF Laboratories

# OOSAFE® IVF LABORATORY DISINFECTANTS

## OOSAFE® DISINFECTANT FOR CO<sub>2</sub> INCUBATORS AND LAMINAR FLOW HOODS



### Excellent material compatibility

- CO<sub>2</sub> incubators
- Laminar flow hoods
- Equipments
- Ultrasound probes
- Hospital furnitures
- Hi-Macs surfaces
- Glass, plastic, acrylic and metal

### Health and safety

- MEA and HSSA tested
- Transparent
- Odorless
- Alcohol-free
- No VOC release

### High quality

- Ready-to-use solution
- Anticorrosive
- Security seal
- Valid for 3 years, even after opening
- Approved by scientific studies
- Long term disinfection/low evaporating rate
- Fast disinfection effect

### Application:

**Spray:** Apply undiluted, spray and distribute evenly on the surface with a sterile cloth. Wait 15 min for the disinfection effect.

**Wipes:** Wipe the surface evenly and wait 15 min for the disinfection effect.

### BACTERICIDAL

Pseudomonas aeruginosa

15 min

Staphylococcus aureus

Enterococcus hirae

### FUNGICIDAL

15 min

Candida albicans

### MYCOBACTERICIDAL AND TUBERCULOCIDAL

15 min

Mycobacterium terrae

Mycobacterium avium

### VIRUCIDAL

15 min

Vaccinia virus\*

Poliovirus

5 min

Adenovirus

TIME UNTIL FULL EFFECT

MICROORGANISM



Product code:  
**OODW-70**

Oosafe® Disinfection Wipes,  
70 wipes per bucket.  
High Quality Wipes:  
260x300 mm, 40 g/m<sup>3</sup>, 100% PP



Product code:  
**OODIH-1000**

Oosafe® Disinfectant for CO<sub>2</sub>  
Incubators and Laminar Flow Hoods,  
1 Liter with spray



Product code:  
**OODIH-5000**

Oosafe® Disinfectant for CO<sub>2</sub>  
Incubators and Laminar Flow Hoods,  
5 Liters refill container



Product code:  
**OODIH-0050**

Oosafe® Disinfectant for CO<sub>2</sub>  
Incubators and Laminar Flow Hoods,  
50 ml in spray bottle

\*Effectiveness against all enveloped viruses including HIV, Hepatitis B and C, Herpes simplex, Influenza A, Ebola

# OOSAFE® IVF LABORATORY DISINFECTANTS

OOSAFE® CLASS IIA MEDICAL DEVICE DISINFECTANT, CE MARKED IN ACCORDANCE WITH EUROPEAN MEDICAL DEVICE DIRECTIVE 93/42/EEC



Manufactured according to the medical device quality management system ISO 13485

## Excellent material compatibility

- CO<sub>2</sub> incubators
- Laminar flow hoods
- Equipment
- Ultrasound probes
- Hospital furniture
- Hi-Macs surfaces
- Glass, Plastic, Acrylic and Metal

## Health and safety

- MEA and HSSA tested
- Transparent
- Odorless
- Alcohol-free
- No VOC release

## High quality

- Ready-to-use solution
- Anticorrosive
- Security sealed
- Valid for 3 years, even after opening
- Approved by scientific studies
- Long term disinfection/low evaporating rate
- Fast disinfection effect

## Application:

**Spray:** Apply undiluted, spray and distribute evenly on the surface with a sterile cloth. Wait 15 min for the disinfection effect.

**Wipes:** Wipe the surface evenly and wait 15 min for the disinfection effect.

### BACTERICIDAL

Pseudomonas aeruginosa

15 min Staphylococcus aureus

Enterococcus hirae

### FUNGICIDAL

15 min Candida albicans

### MYCOBACTERICIDAL AND TUBERCULOCIDAL

15 min Mycobacterium terrae

Mycobacterium avium

### VIRUCIDAL

15 min Vaccinia virus\*

Poliovirus

5 min Adenovirus

TIME UNTIL FULL EFFECT

MICROORGANISM



Product code:  
**OOMDW-120**  
Oosafe® Disinfectant Wipes,  
70 wipes per bucket.  
High Quality Wipes:  
260x300 mm, 40 g/m<sup>2</sup>, 100% PP



Product code:  
**OOMD-01000**  
Oosafe® Disinfectant,  
1 Liter with spray



Product code:  
**OOMD-05000**  
Oosafe® Disinfectant,  
5 Liters refill container



Product code:  
**OOMD-00050**  
Oosafe® Disinfectant,  
50 ml in spray bottle

\*Effectiveness against all enveloped viruses including HIV, Hepatitis B and C, Herpes simplex, Influenza A, Ebola

# OOSAFE® IVF LABORATORY DISINFECTANTS

## OOSAFE® HAND DISINFECTANT



### Health and safety

- MEA and HSSA tested
- Transparent
- Odorless
- Alcohol-free
- No VOC release
- Keeps your skin soft

### High quality

- Ready-to-use solution
- Security sealed
- Valid for 3 years, even after opening
- Long term disinfection/low evaporating rate
- fast disinfecton effect

### Application:

Rub hands with 2 pumps (3 ml) disinfectant for at least 30 seconds until hands are completely dry.

TIME UNTIL FULL EFFECT	MICROORGANISM	
	BACTERICIDAL	FUNGICIDAL
30 sec	Pseudomonal aeruginosa	
	Staphylococcus aureus	
	Enterococcus hirae	
30 sec	Escherichia coli	
		Candida albicans
30 sec		Vaccinia virus*



Product code:  
**OODH-0500**  
Oosafe® Hand Disinfectant,  
500 ml with pump



Product code:  
**OODH-5000**  
Oosafe® Hand Disinfectant,  
5 Liters refill container



Product code:  
**OODH-0050**  
Oosafe® Hand Disinfectant,  
50 ml in spray bottle

\*Effectiveness against all enveloped viruses including HIV, Hepatitis B and C, Herpes simplex, Influenza A, Ebola

# OOSAFE® IVF LABORATORY DISINFECTANTS

## OOSAFE® LABORATORY SURFACE DISINFECTANT



### Health and safety

- MEA and HSSA tested
- Clear - slightly yellow
- Alcohol-free
- Very low level of VOC release


### High quality

- Security sealed
- Valid for 3 years, even after opening
- Long term disinfection

### Excellent material compatibility

- Laboratory surfaces and floor
- Glass, Plastic, Tiles, Rubber

**Application:** Dilute according to the table below with water and after application wait 30 min for the required disinfection effect. Please rinse applied area after application.

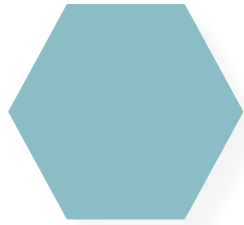

Product code:  
**OODSF-02000**  
Oosafe® Laboratory  
Surface Disinfectant,  
2 Liters with dose cup



Product code:  
**OODSF-10000**  
Oosafe® Laboratory  
Surface Disinfectant,  
10 Liters refill container



Product code:  
**OODSF-00050**  
Oosafe® Laboratory  
Surface Disinfectant,  
50 ml in bottle

MICROORGANISM		DILUTION
<b>BACTERICIDAL</b>		
15 min	Pseudomonas aeruginosa	1:50
	Staphylococcus aureus	1:50
	Enterococcus hirae	1:50
<b>FUNGICIDAL</b>		
15 min	Candida albicans	1:50
<b>MYCOBACTERICIDAL AND TUBERCULOCIDAL</b>		
30 min	Mycobacterium terrae	1:25
	Mycobacterium avium	1:25
<b>VIRUCIDAL</b>		
15 min	Vaccinia virus*	1:50
30 min	Poliovirus type 1	1:25
10 min	Adenovirus type 5	1:25

\*Effectiveness against all enveloped viruses including HIV, Hepatitis B and C, Herpes simplex, Influenza A, Ebola



THE DANISH CHILDHOOD  
CANCER ORGANISATION

SparMED ApS - Ryttermarken 2, 3520 Farum, Denmark - Phone (+45) 39 40 25 03 - [www.sparmed.dk](http://www.sparmed.dk)



 [facebook.com/sparMED](https://facebook.com/sparMED)

 [twitter.com/sparMED](https://twitter.com/sparMED)

 [linkedin.com/company/sparmed-aps](https://linkedin.com/company/sparmed-aps)